

# Fully automated artifact reduction method for time-resolved cone-beam CT angiography

Chengzhu Zhang<sup>1</sup>, Xin Tie<sup>1</sup>, Yijing Wu<sup>1</sup>,  
John W Garrett<sup>1</sup>, Guang-Hong Chen<sup>1,2</sup>

1. Department of Medical Physics, University of Wisconsin, Madison, WI
2. Department of Radiology, University of Wisconsin, Madison, WI



DEPARTMENTS OF  
**Medical Physics & Radiology**  
UNIVERSITY OF WISCONSIN SCHOOL OF MEDICINE AND PUBLIC HEALTH



# Acknowledgement


- This work is supported by NIBIB of NIH under award U01 EB021183



# Motivation: one-stop-shop (OSS) protocol

- Large Vessel Occlusion (LVO) stroke patients benefit from multi-modal imaging.
- Despite the benefits of performing this multi-modal imaging, the time needed to perform this imaging process results in significant delays to patient treatment<sup>1</sup>.


- The desire to save time has motivated the one-stop-shop imaging paradigm for acute stroke imaging<sup>2</sup>.

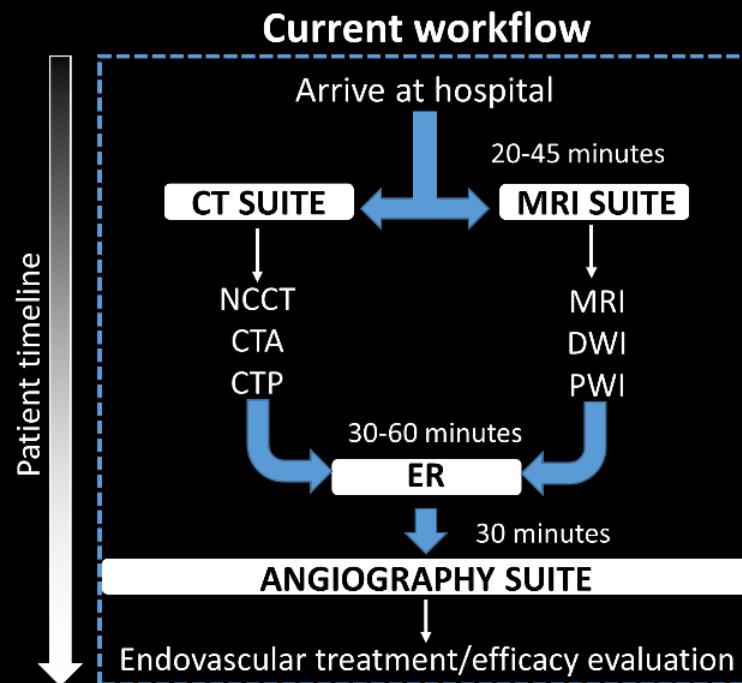


potentially save at least  
one hour per patient!



# Introduction: one-stop-shop (OSS) protocol

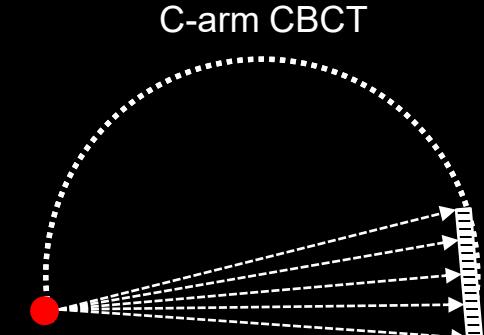




eSMART-  
RECON<sup>1</sup>

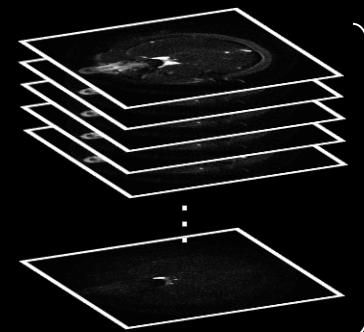

1. Li et al, IEEE TMI (2019)


256x256x200  
x 40 frames

IV  
Injection

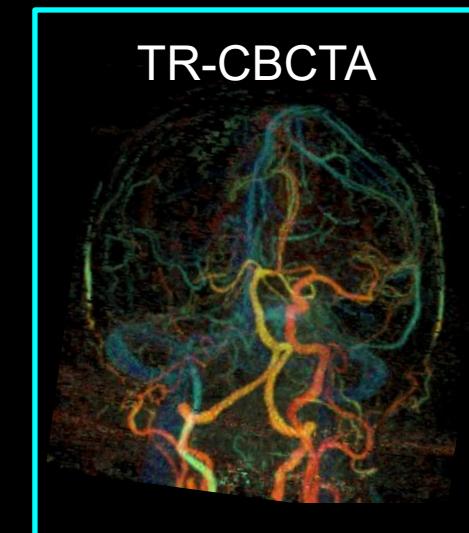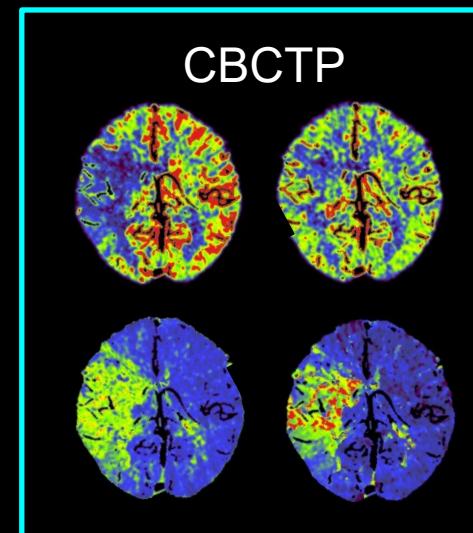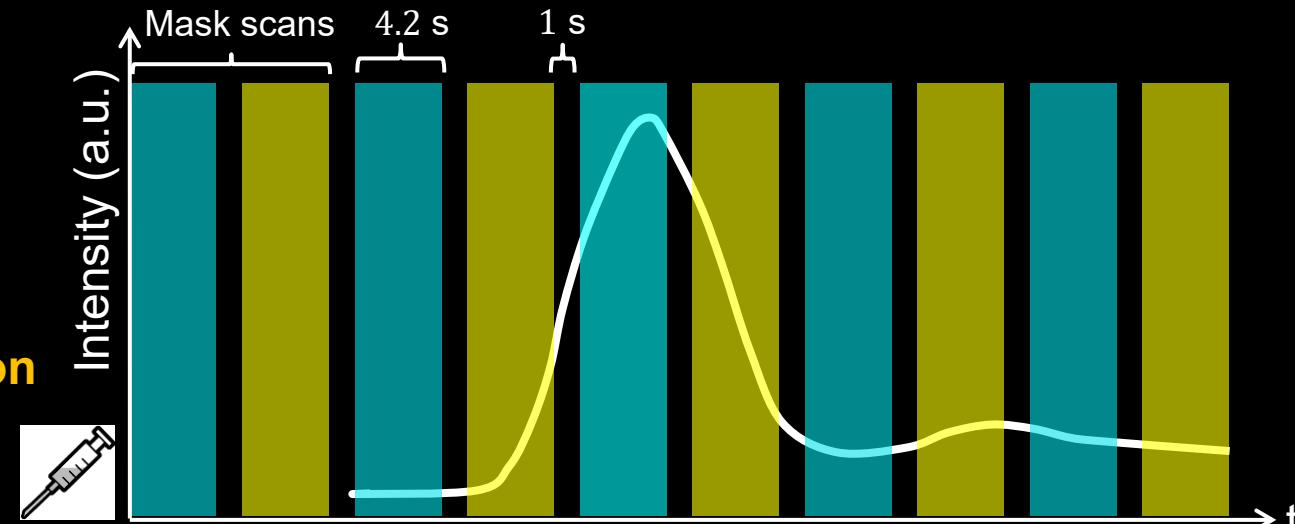

# Introduction: the benefit of TR-CBCTA


- Besides collateral analysis, color-coded time-of-arrival maps<sup>1</sup> can be derived from TR-CBCTA and provide complementary information of various vasculature phases.

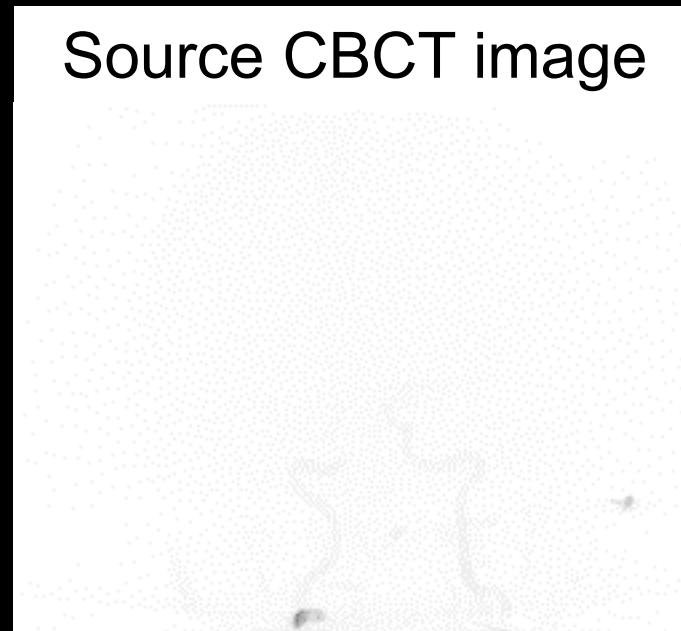

Time @ intensity reaches 30% peak



Time-of-arrival map

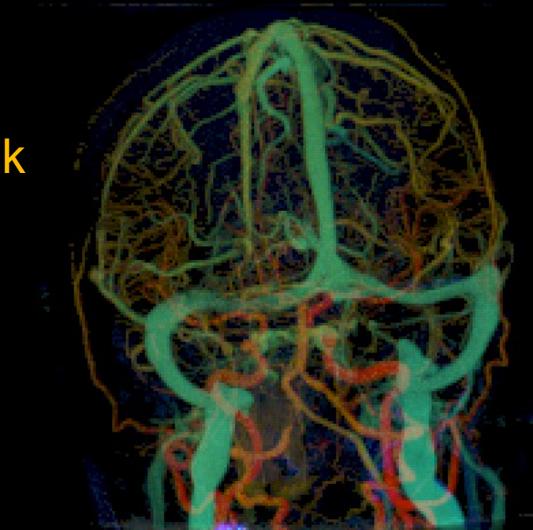

Early arrival




Late arrival




# The technical limitation

- Time-resolved cone-beam CT angiography (TR-CBCTA) acquired in the angiography suite often presents **mis-registration** and **residual bone artifacts** due to patient motion and gantry instability.

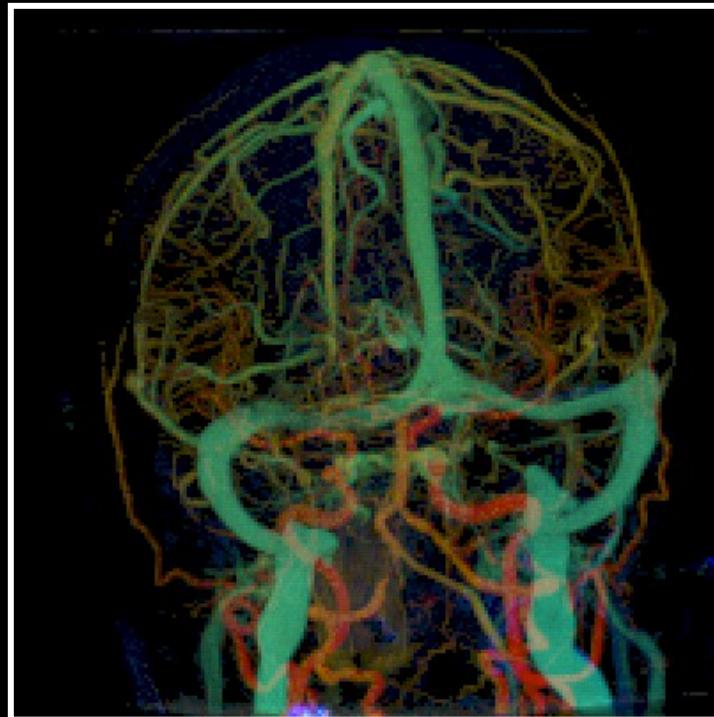


**Ideally (very few)**

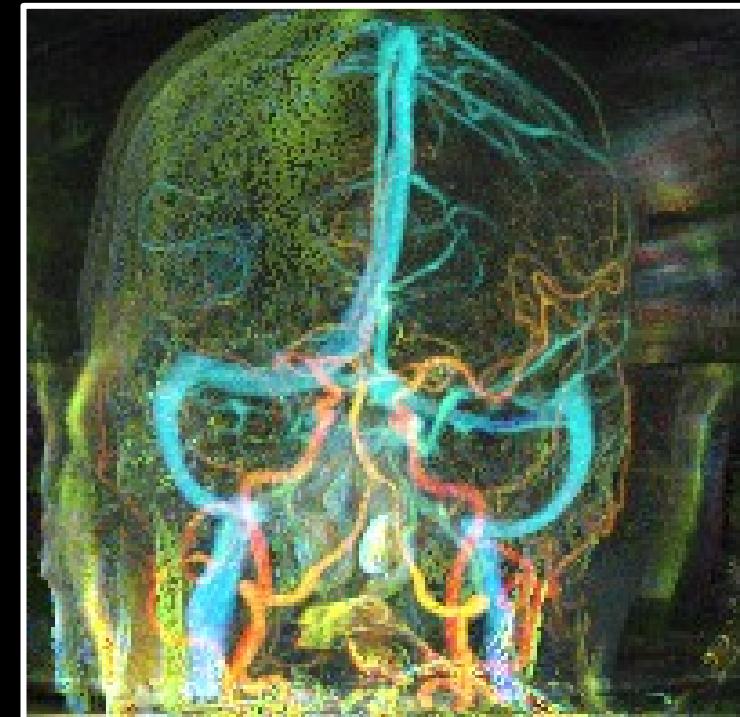

**In reality (majority)**



# Conventional supervised deep learning method

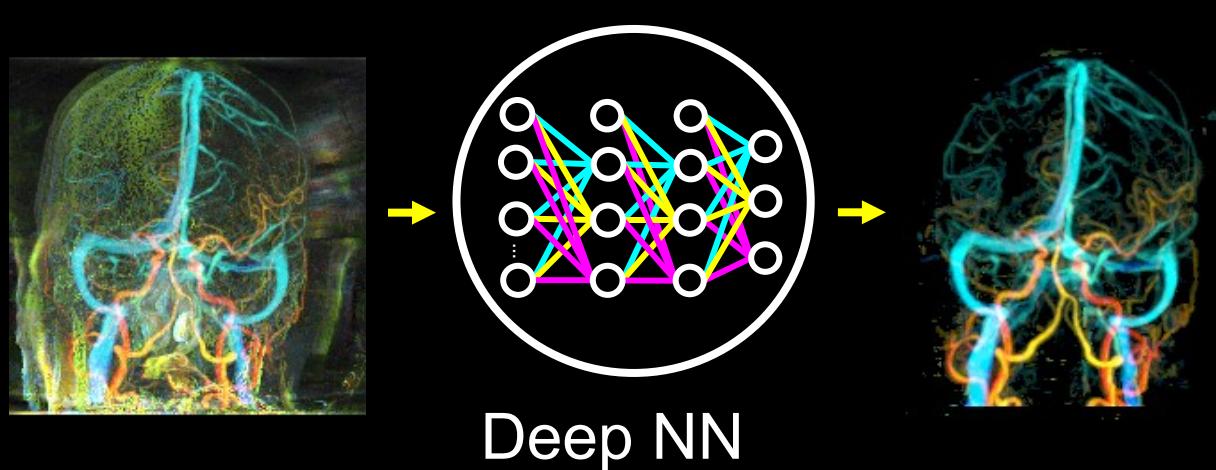

- Implementation details<sup>1</sup>
  - Network architecture: U-Net
  - Input: 256 x 256 (axial image)  
x 40 frames, TR-CBCTA  
via eSMART-RECON<sup>2</sup>
  - Training target: 256x256 hand-labeled vasculature map
  - Output: 256x256 Probability map
  - Training data: 50 patient cases,  
5-fold cross-validation




# Laborious hand annotation process

- Vasculature is hand-labeled in Temporal-MIP images for each patient.
  - >200k labeled voxels/patient
- Need special medical domain knowledge.

Temporal-MIP CTA

Labels

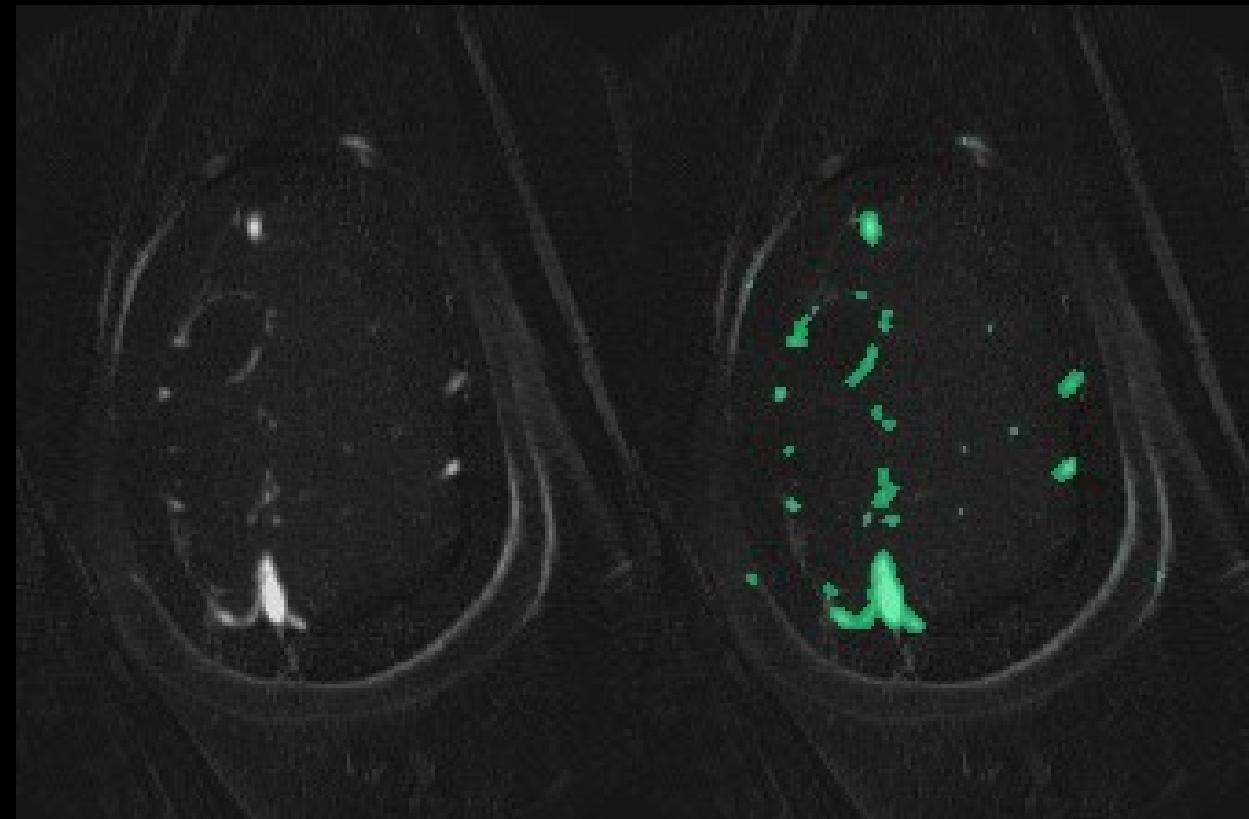




# Objective

- To develop a **fully-automated** machine learning method to effectively reduce artifacts in TR-CBCTA in angiography suites.

# Method: hybrid artifact reduction framework

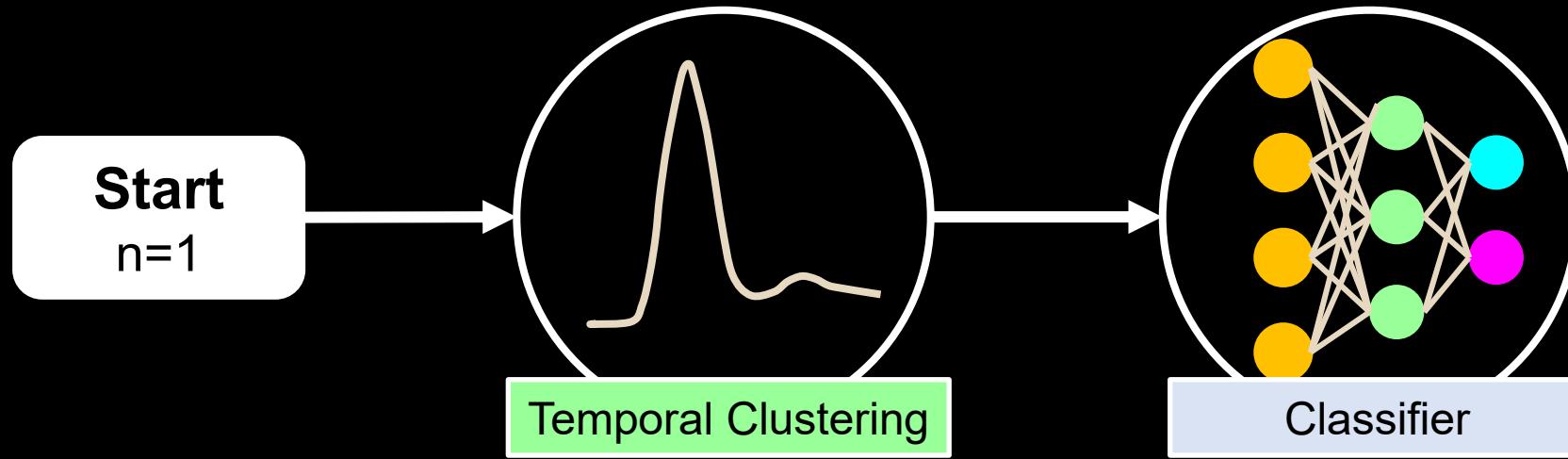


Input: SMART-RECON  
Denoised CTA<sup>1,2</sup>  
Volume dim: 256x256x200  
Time points: 40

Method: K-Means clustering  
Input: 40-D temporal curves  
Output: 4 clusters  
Format: volumetric mask

Method: Supervised Learning  
Input: 256x256x200 clustering masks  
Target: manual labeling, 'vessel', 'artifact'  
Output: prediction score (threshold=0.5)




# Method: artifact reduction framework

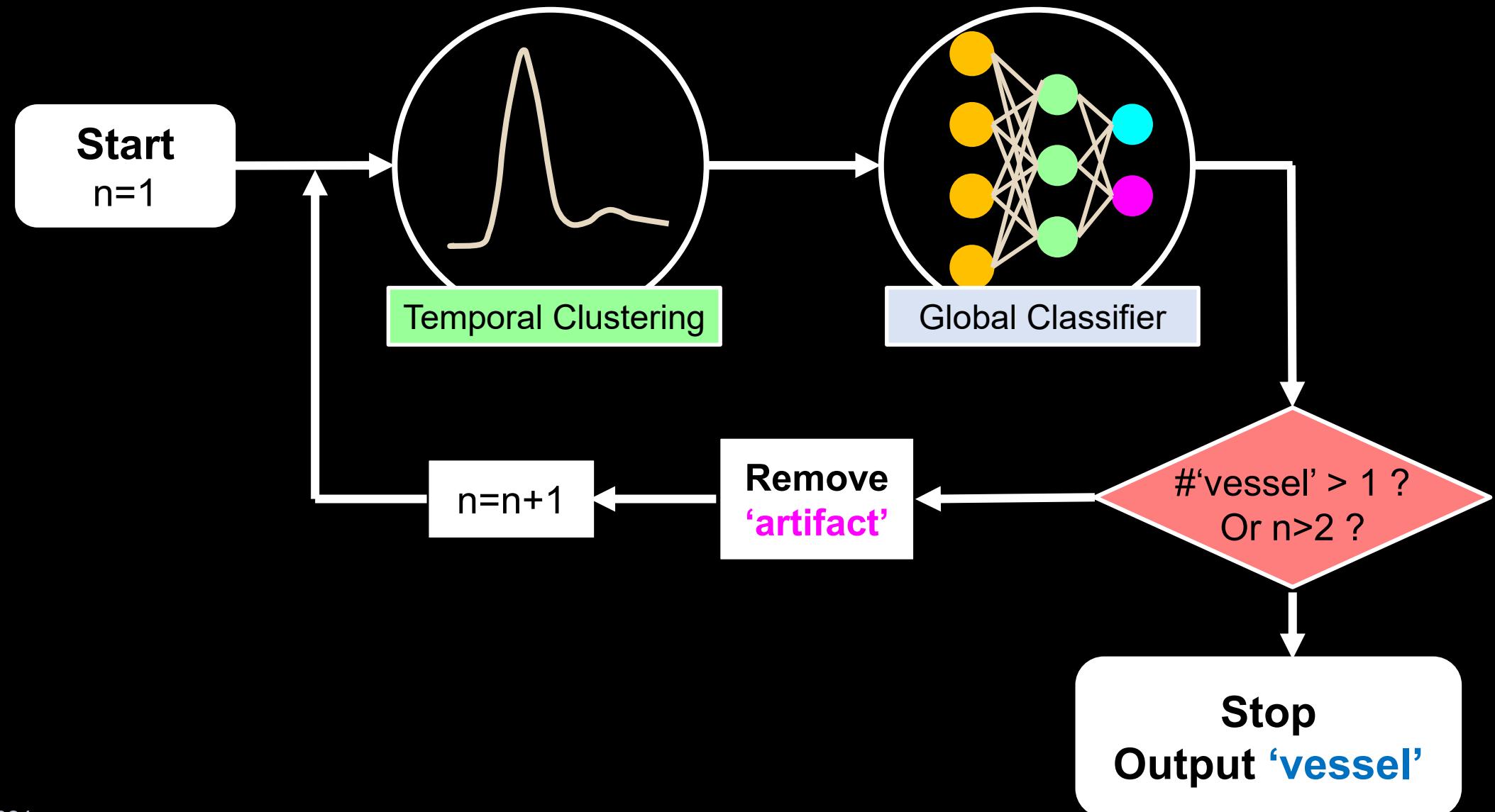




# An example of working pipeline

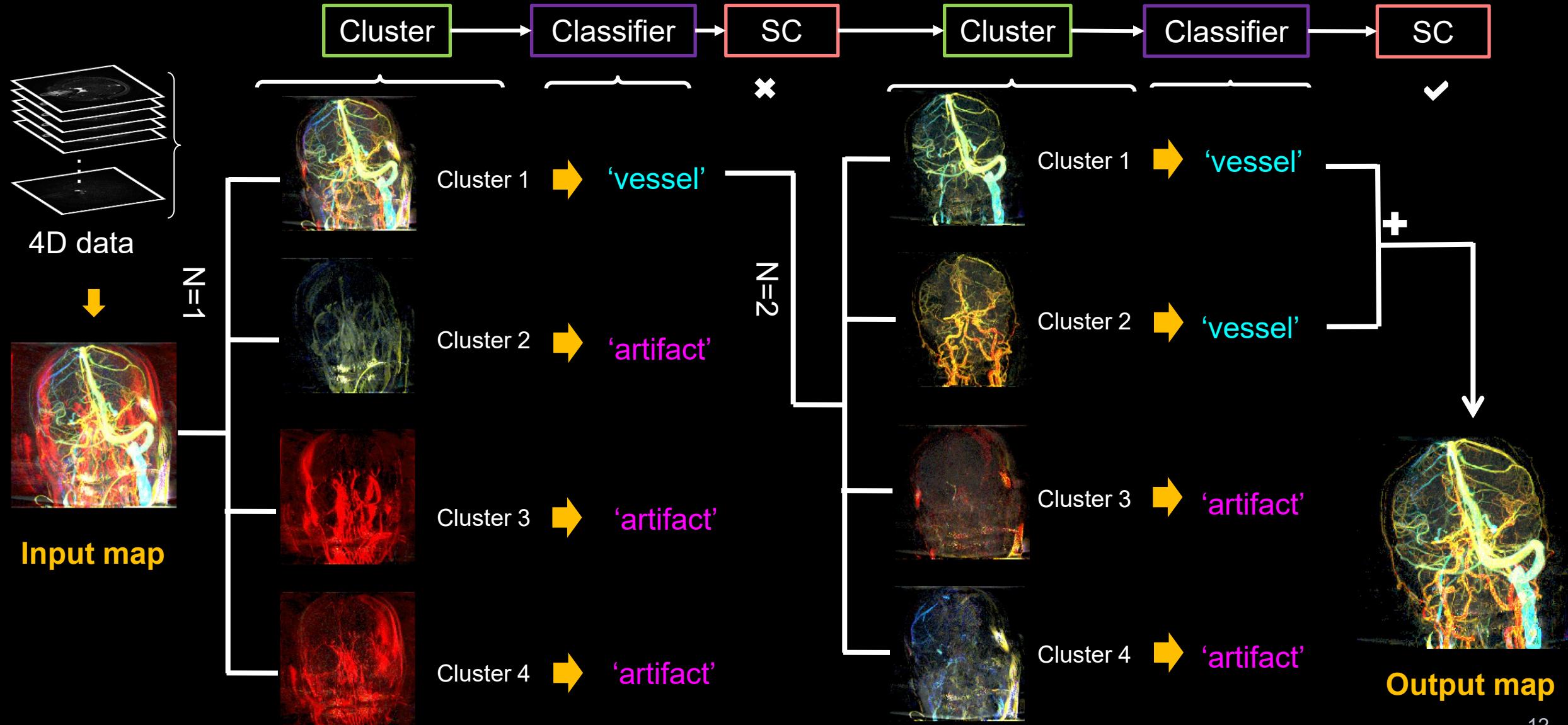




# Methods: study schema

- IRB approved and HIPAA compliant study

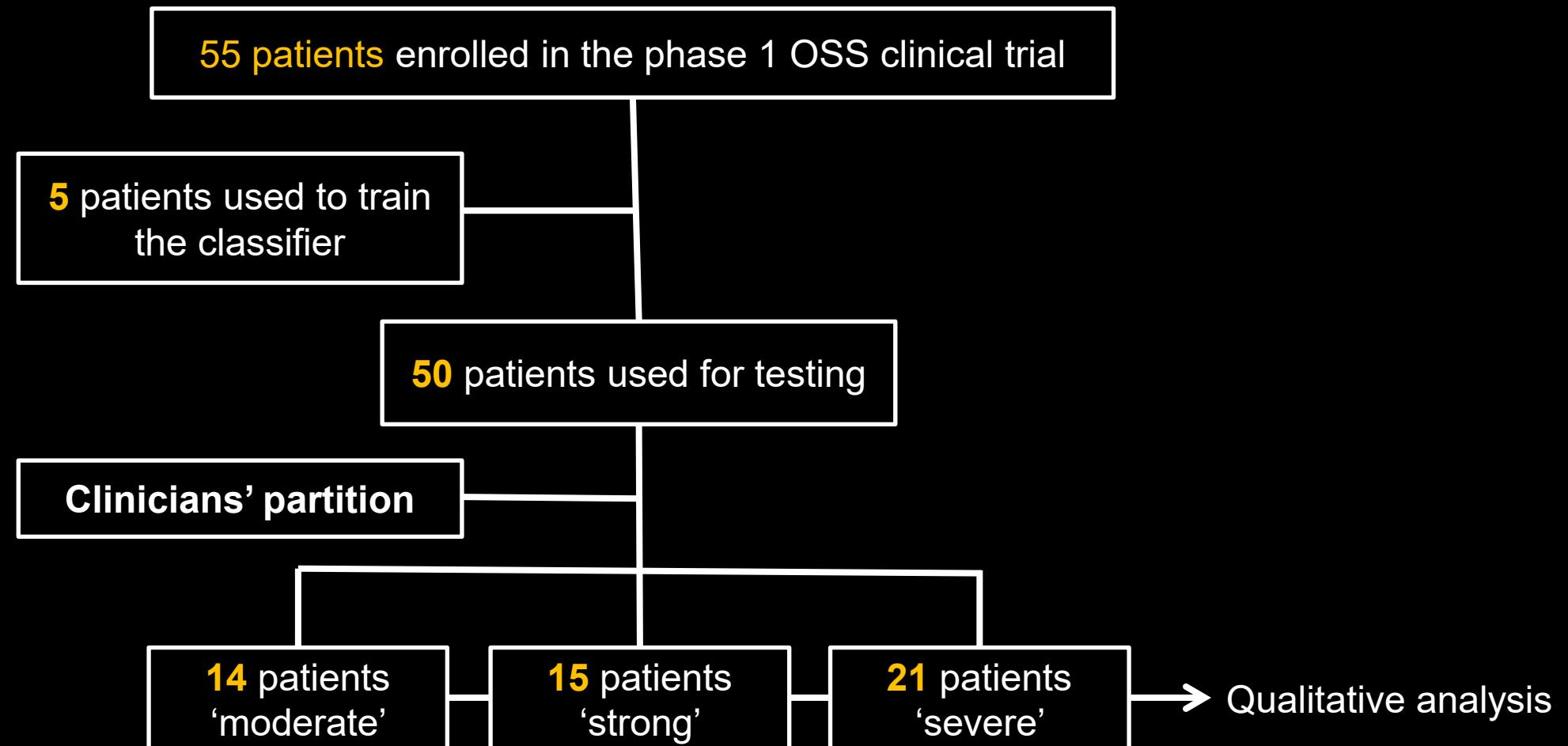


# Results I: a case with strong artifacts

Original

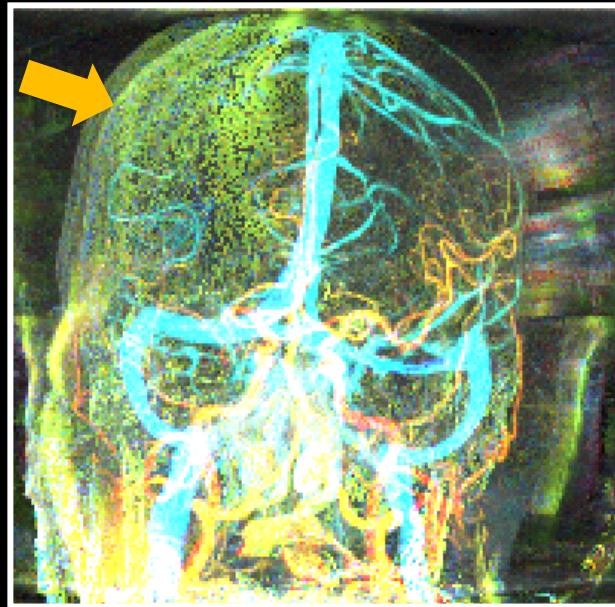


Proposed

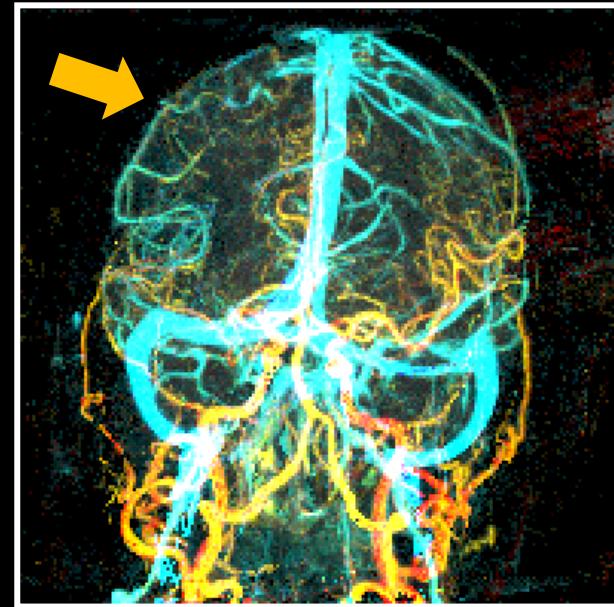


# Results II: a case with strong artifacts


Original

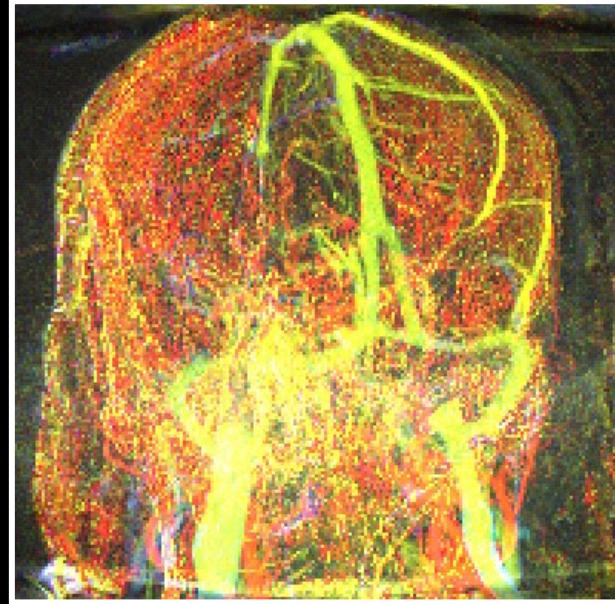


Proposed

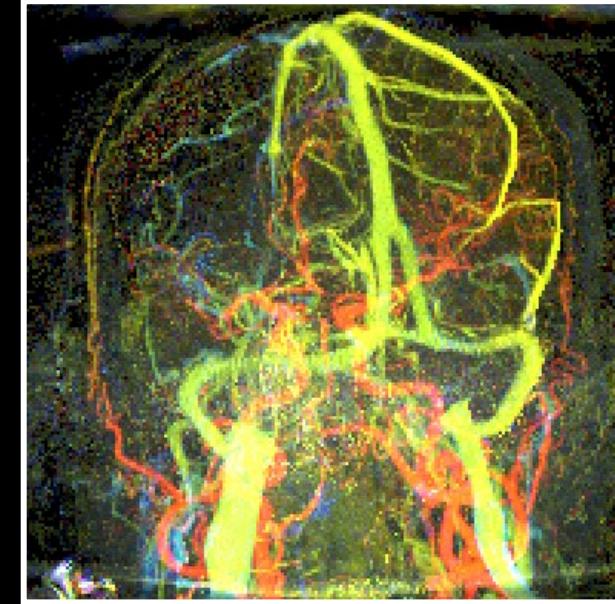

'Artifact'

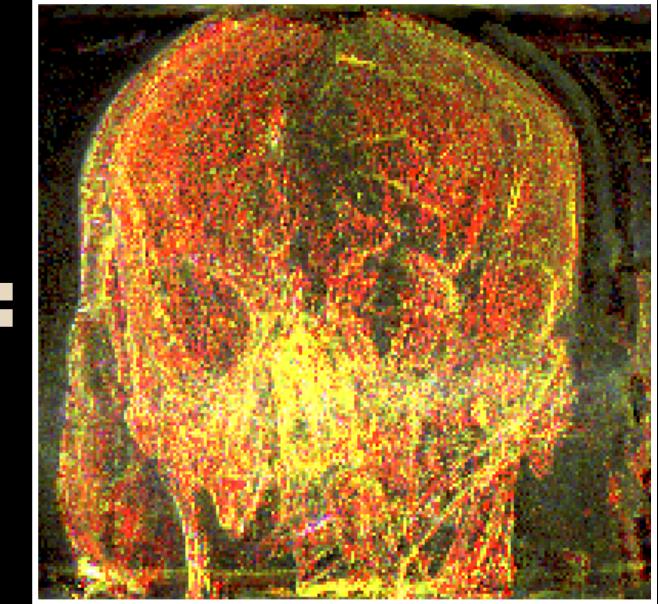


# Results III: a case with severe artifacts

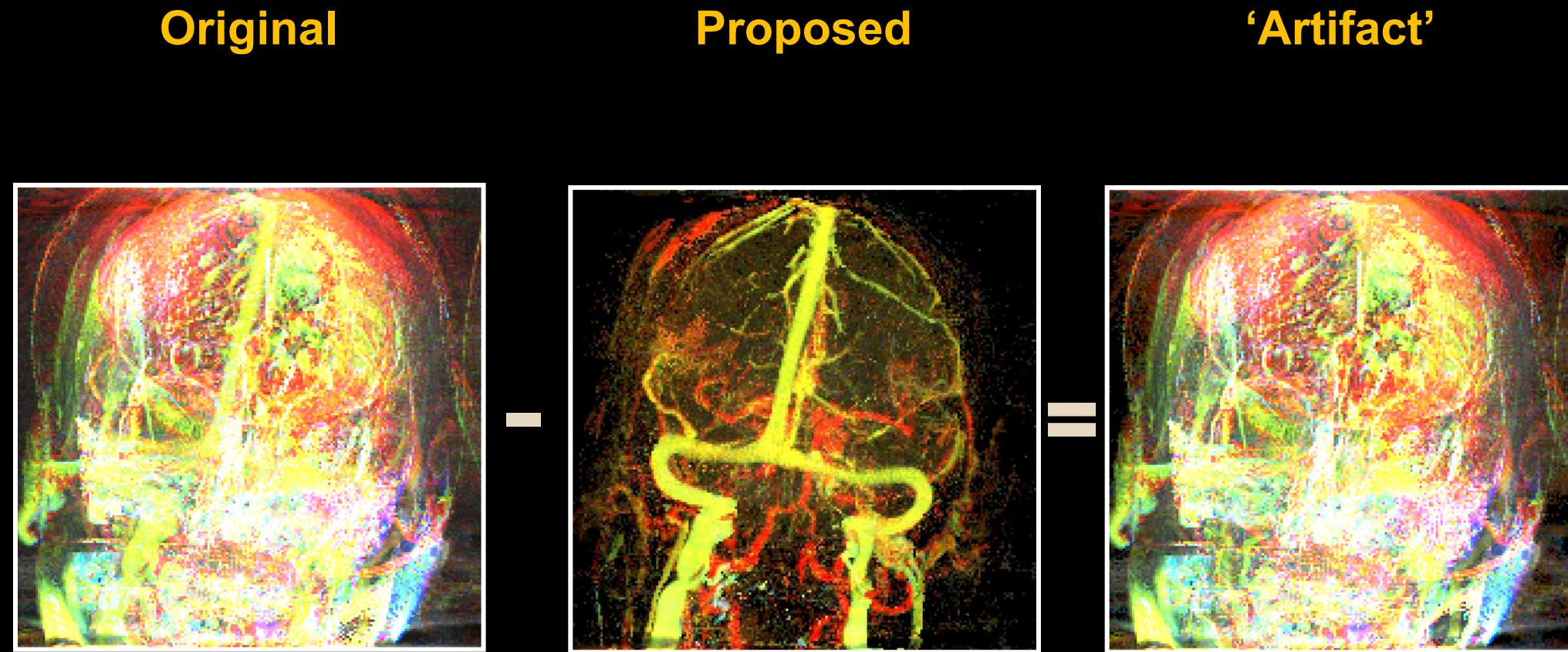


# Robustness analysis: moderate artifacts

DeepNN<sup>1</sup>

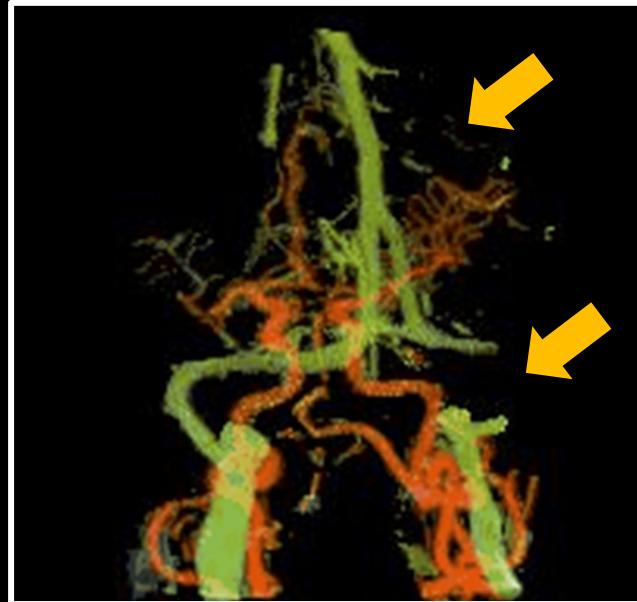


Proposed

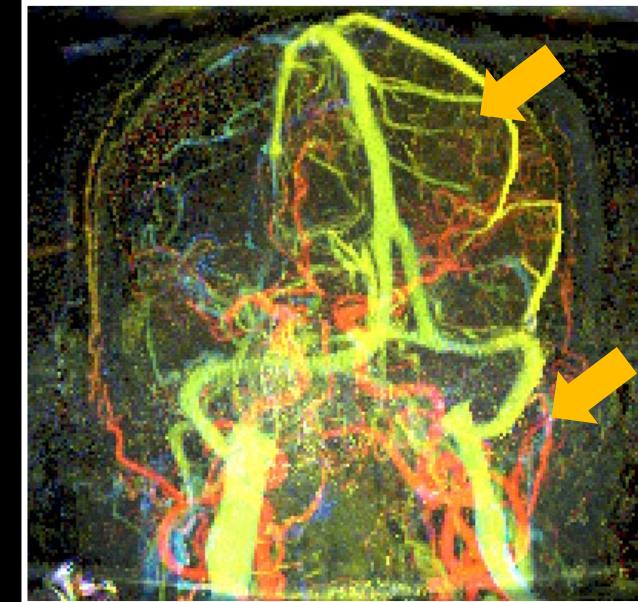


1. Juan et al; RSNA. (2019)

# Robustness analysis: severe artifacts

DeepNN<sup>1</sup>

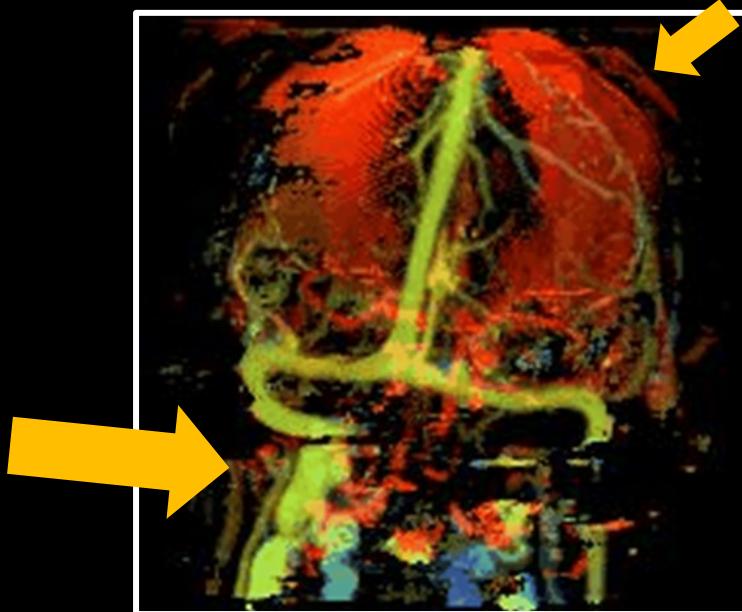


Proposed

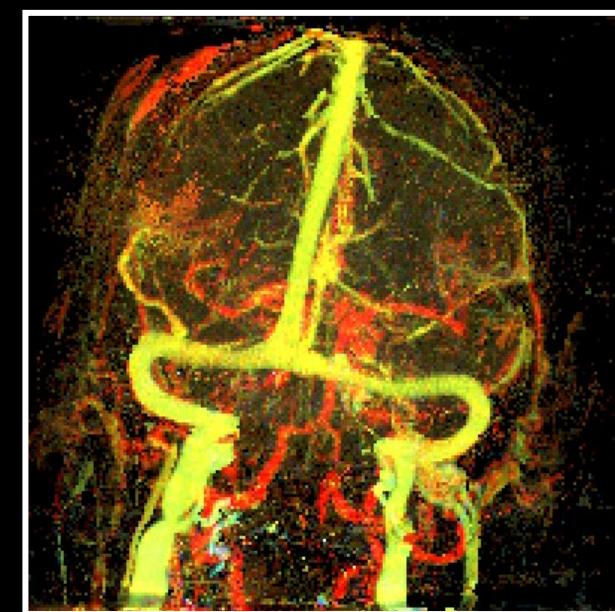


1. Juan et al; RSNA. (2019)

# Conclusion

- Robust artifact reduction has been achieved for clinical cases in the test cohort regardless of their initial artifact levels.
- After artifact reduction, the visualization of TOA maps clearly shows different contrast enhanced phases of vessels.

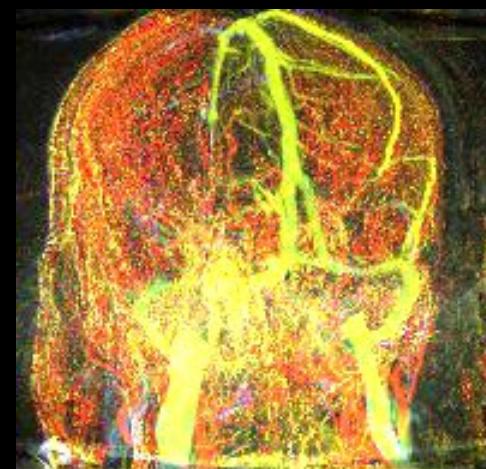

Proposed

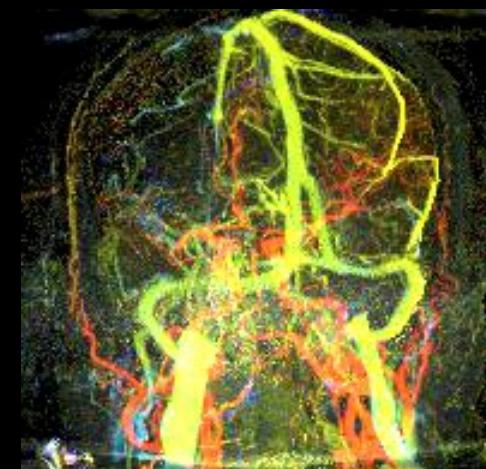



# Acknowledgement

Thank You  
[cchang553@wisc.edu](mailto:cchang553@wisc.edu)

